


External power adapter
Avoid security risks

Free flow clamp
Self-flow prevention

Data interconnection
Display infusion therapy and sign
monitoring data in one screen

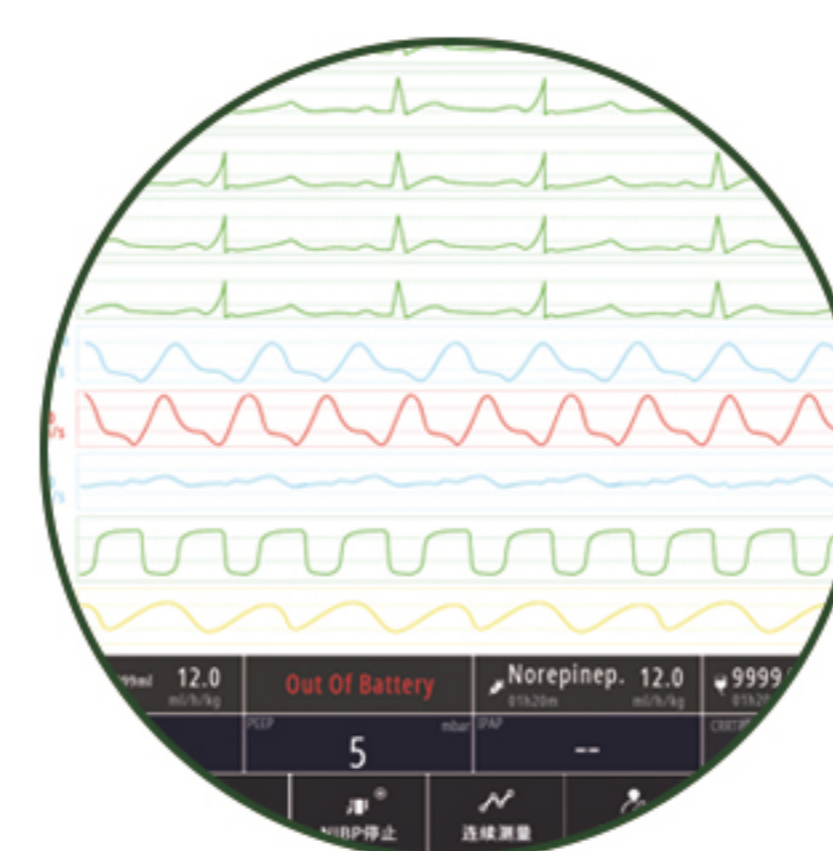

MI27



Safe and error proof



Smart Infusion



Flexible transfer

Data interconnection

Support HL7 protocol
infusion log data is open to share, collect the data in time and
rejection of information islands

Speed up starting the machine

Start up one second, every minute counts

DERS

Drug library support: Built-in drug list and DERS,
alarm the error of drug and dose, ensure the safety of drug use

Safe infusion

Triple dose anti-stall design, ten safety measures for infusion



Free combination of
multiple channels

Infusion workstation

Infusion Pump MI27

PRODUCT SPECIFICATION

Dimensions	216×131×72mm (width×depth×height)
Weight	1.26kg
Screen	MI27: 4-inch HD display
Classification	Class II, Type CF, IP24

BAISC PARAMETERS

Volume infused	0.01-9999.99ml, step by 0.01ml
Flow rate	0.01-2000ml/h, step by 0.01ml/h
Time infused	00h00min~99h59min, step by 1 min
Accuracy	±4.8% ±2% (MDK Infusion set)
Purge / Bolus rate	1ml/h~2000ml/h; step by 1ml/h volume 0.01~9999.99ml adjustable
Infusion Mode	RVT, Dose, Weight, Relay, Gradient, Drip, First Delivery, Drug, Loading Dose, Micro, Sequence, Intermittent, etc
DERS	DERS intelligent medication system
Information Transmission	Support WIFI、HL7、Wired
KVO flow rate	Adaptable or Automatic
Alarm	18 kinds of alarms
Occlusion pressure	3 levels ,Maximum: 120 kPa±10kPa, Minimum: 20 kPa±10kPa
Bubble sensitivity	Minimum detectable≥5μl Bubble, 8-speed adjustable

POWER SUPPLY

Power Adaptor	A.C. input 100V-240V 50/60Hz 1.0A D.C. output 15V 2A max
Battery	10.8V rechargeable Lithium battery
Battery capacity	Support normal operation for 9 hours

MDKMed Medical Technology Co.Ltd.

Tel: 4008-808-392

Email: sales@mdkmed.com

Website: www.mdkmed.com

Company Address: Floor 20, Poly Center, shangcheng District, Hangzhou City, Zhejiang, China.

Registered Company: Zhejiang MDKkingdom Technology Co., Ltd.

MDKMed is the registered trademark of MDKMed Medical Technology Co., Ltd. Specifications are subject to change without prior notice.
©2021 MDKMed Medical Technology Co., Ltd. Copyright Reserved.



Learn more