SHOCKWAVE MACHINE

USER MANUAL



Model No :SW14



Model No :SW15

Warning

Dear user:

Thanks for using our company's products, in order to fully make use of the products, we sincerely suggest follow:

Please read the instruction in detail and keep it well for further reading and consulting.

◆◆Please follow the instruction to fit and operate the machine correctly.

◆◆Please do not remove and change any accessories of the machine.

◆◆It forbids others to open and remove the mainframe except the authorized people in the company.

◆◆It forbids others to change and maintain the machine except the professional people.

◆◆Please contact us if there is any trouble during operation.

General Introduction

The shockwave system uses the ballistic principle of shockwave generation: A pressure wave is formed via a projectile by using accelerated compressed air. The compressed air is generated by an electronically-controlled ballistic-pressure compressor. Using elastic impact, the kinetic energy of the projectile is transferred into the probe of the applicator and then into the client's body. Consequently, during the treatment, the end of the applicator must be in direct-contact with the skin and subcutaneous tissue.

Shockwave is aimed at the affected areas that are the source of chronic pain. The influence of the shockwave causes to the dissolution of calcium deposits and leads to better vascularization. The after-effects are relief from the pain.

Shockwave has the following effects:

Cellular: Increase in cell membrane transmittance by improving ionic channels activity, stimulation of cell division, stimulation of cellular cytokines production.

Reproduction of vessels in the area of tendons and muscles: Improvement of blood circulation, increase in concentration of growth factor beta 1, chemotactic and mitogenic effect on osteoblasts.

Effect on nitrogen oxide system: Bone healing and remodeling.

Improvement of macro-circulation and metabolism.

Dissolution of calcified fibroblasts.

Supports the production of collagen.

Reduction in tissue tension.

Analgesic effect.

Sw14 Model Interface introduction:



a

1. **Power+**
2. **Power-**
3. **Forearm mode, frequency 4HZ, output level 35MJ**
4. **Upper arm mode, frequency 6HZ, output level 55MJ**
5. **Thigh mode, frequency 6HZ, output level 55MJ**
6. **Leg mode, frequency 5HZ, output level 30MJ**
7. **Buttock mode, frequency 7HZ, output level 105MJ**
8. **Back mode, frequency 16HZ, output level 50MJ**
9. **Shockwave woke head jack**
10. **Frequency+**
11. **Frequency-**
12. **Abdomen mode, frequency 10HZ, output level 25MJ**
13. **Shoulder mode, frequency 8HZ, output level 50MJ**
14. **Twist mode, frequency 12HZ, output level 75MJ**
15. **Scapular mode, frequency 3HZ, output level 105MJ**
16. **Low frequency mode, frequency 1HZ, output level 105MJ**
17. **general-purpose mode**
18. **Power interface&Power switch(At the bottom of the suitcase)**

Tips: The working times is set in 2500 shots, after that the machine will stop to output energy. And then press the begin button the machine will start work again.

The machine is acquiesce in general-purpose mode

（Only in general-purpose mode，the frequency level and power level can be adjusted. Below is the details of the general-purpose mode）

Sw15 Model Interface introduction:



1. **Power+**
2. **Power-**
3. **Shoulder mode, frequency 8HZ, output level 50MJ**
4. **Upper arm mode, frequency 6HZ, output level 55MJ**
5. **Forearm mode, frequency 4HZ, output level 35MJ**
6. **Thigh mode, frequency 6HZ, output level 55MJ**
7. **Leg mode, frequency 5HZ, output level 30MJ**
8. **general-purpose mode**
9. **Shockwave woke head jack**
10. **Frequency+**
11. **Frequency-**
12. **Scapular mode, frequency 3HZ, output level 105MJ**
13. **Twist mode, frequency 12HZ, output level 75MJ**
14. **Abdomen mode, frequency 10HZ, output level 25MJ**
15. **Back mode, frequency 16HZ, output level 50MJ**
16. **Buttock mode, frequency 7HZ, output level 105MJ**
17. **Low frequency mode, frequency 1HZ, output level 105MJ**
18. **Power interface&Power switch(At the bottom of the suitcase)**

Tips: The working times is set in 2500 shots, after that the machine will stop to output energy. And then press the begin button the machine will start work again.

The machine is acquiesce in general-purpose mode

（Only in general-purpose mode，the frequency level and power level can be adjusted. Below is the details of the general-purpose mode）

General parameter introduction

Output frequency gears in general mode:

1=1HZ 2 =2HZ, 3 =3HZ,4=4HZ,5=5HZ,6=6HZ, 7=7HZ,

8 = 8HZ, 10 = 10HZ, 11 = 11HZ, 12 = 12HZ, L3 = 13HZ, 14 = 14HZ,

15 =15HZ, 16 =16HZ

Output level Grade of the general-purpose mode: 1-12

Output level: When Frequency in 1Hz-12HZ, 1=25MJ, 2=30MJ, 3=35MJ, 4=40MJ, 5=45MJ,

6=50MJ, 7=55MJ, 8=65MJ, 9=75MJ, 10=85MJ, 11=95MJ, 12=105MJ

When Frequency in 13Hz-14HZ, 1=25MJ, 2=30MJ, 3=35MJ, 4=40MJ, 5=45MJ, 6=50MJ,

7=55MJ, 8=60MJ, 9=65MJ, 10=70MJ, 11=75MJ, 12=80MJ


When Frequency in 15Hz-16HZ, 1=25MJ, 2=28MJ, 3=31MJ, 4=35MJ, 5=38MJ, 6=41MJ,

7=45MJ, 8=48MJ, 9=50MJ, 10=52MJ, 11=55MJ, 12=58MJ

Suggested value

Application Treatment Recommendation					
Area	**Frequency/P**	**shocks**	**Sessions/time**	**total shocks**	**Interval**
Elbow	3HZ/P1	2000-2500	3-5	6000-12500	5-7days
Hand/Wrist	1HZ/P1	2000	3-5	6000-10000	5-7days
Pubic Symphsis	3HZ/P3	2000-2500	5	6000-12500	every day/ 1 day
Knee	3HZ/P1	2000	3-5	6000-10000	5-7days
Foot/Ankle	3HZ/P1	1000	3-5	3000-5000	1days
Shoulder	6HZ/P5	2000-2500	3-5	6000-12500	5-7days
Lumbar Region	8HZ/P7	2000-2500	3-5	6000-12500	5-7days
Fat accumulate	10HZ/P7	2000-2500	5	6000-12500	5-7days
Hip	4HZ/P10	2000-2500	5	6000-12500	5-7days

Handles

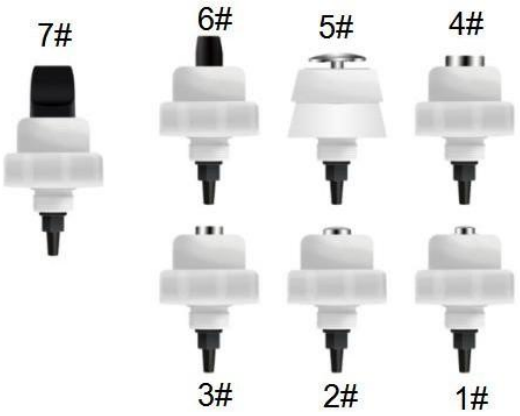



Tip:no.6 and 7 for ED treatment,1.2.3.4for body arm leg treatment 5 for prostate

The blue cover is used with the 5 tips.



You can buy extra the gray hard tips,not included inside the machine.

#1 is for ED treatment
#2 is for elbow treatment
#3 is use together with the tip 5

For ED Treatment advice


1. The frequency of ED treatment is 1-10HZ. The recommended frequency is 3-8HZ.
2. The energy of ED **advice frequency:2-3**
P1-P4
and the energy is adjusted according to the size of the customer.
3. Treatment time: 2 times a week for 3 weeks, 3 weeks after the break, the second treatment starts at the 7th week
4. Treatment location. Along the 5 areas of the spongy tissue, each area was played about 300 times. That is, 1500 treatments per time, about 20 minutes
5. You can play several times at the far end of the penis and play it a few times at the proximal end of the root.

Most men suffering from erectile dysfunction have vascular problems affecting the vessels that supply blood to the cavernous bodies of the penis, resulting in a decreased ability to develop and maintain an erection. The end result is an inadequate blood-flow followed by premature draining of an already reduced amount of blood from the cavernous bodies, meaning an erection is not achieved. Shockwave Therapy for ED of this type can be a highly effective treatment.

ED Therapy Protocol

1. Target to spongy tissue in 5 different areas.
2. 300 impacts per area, total 1500 impacts in one treatment.
3. Twice a week , treat consistently for 3 weeks, then stop for 3 weeks. And start next course .
4. Should be more impacts at the top end of penis and less impacts at the bottom of penis.
5. Operation guide:


(1) **Choose the 6# and 7# exclusive ED treatment work head**, Along the 5 areas of the spongy tissue, 300 impacts per area, total 1500 impacts in one round, about 20 minutes.

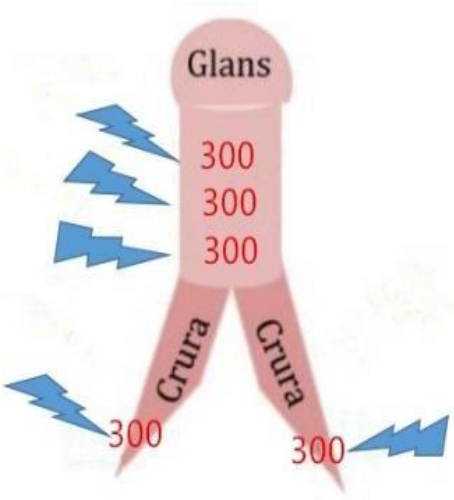








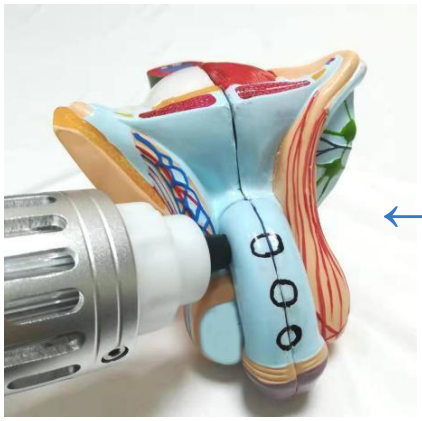

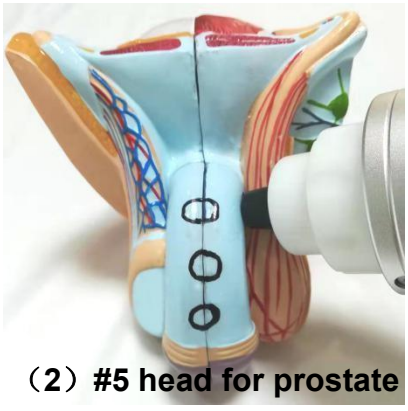

（2）**#5 head for prostate**







Shockwave for Cellulite Reduction

Shockwave, or Acoustic Radial Wave Therapy works by stimulating the fat breakdown, collagen synthesis, and lymphatic drainage.

Our shockwave is delivered by placing a probe up against the areas of cellulite, as it passes over the area it emits radial waves through the skin.



Shock wave for Physiotherapy

Shock wave is an acoustic wave which carries high energy to painful areas and soft tissues with subacute, subchronic and chronic conditions. This energy promotes healing, regenerating and reparative processes. It's a unique, non-invasive solution for pain associated with the muscul oskeletal system.

Three to four treatments needed at weekly intervals. A therapy session only takes approximately 10 minutes.

Warranty description

According to the company's regulations, the machine warranty period is one year, and the warranty period for the handles and accessories is 3 months.

1. During the warranty period, if it is a fault caused by the quality of the product itself, please contact the after-sales department and present this warranty card.
2. Do not accept machines that have failed due to unauthorized modifications.
3. Once the warranty card has been altered, the warranty will be invalid immediately.
4. Please keep the warranty card as the warranty certificate of this machine.

The following conditions are not free of repair.

1. No warranty certificate
2. Failure caused by do not operating the machine in accordance with the instructions or video.
3. Damage caused by disassembling the machine without the permission of the company.
4. Failure, scratching or damage due to movement or falling.
5. The user has improper storage and damage caused by improper use.

Warranty card

Product Model: _____ Product Name: _____

Purchase Store: _____ Purchase Date: _____

Customer Name: _____ Customer Phone: _____

Purchase Store: _____ Purchase date: _____

Product fault description:_____

_____

_____

_____

_____